

Gel-Pak™ Gel Compresses

Instant Cold Compresses


Therapeutic pack for the use with common aches, pains, swelling and bruising. Perfect to keep in your First Aid kit or in your sport bag.

Characteristics

- Instant cold chemical reaction
- Single use
- Stays cold approximately 30 minutes
- Non-toxic product

Available sizes

CODE	DESCRIPTION	DIMENSIONS
712	COLD THERAPY COMPRESS WITH SOFT TEXTURE (CAN BE APPLIED DIRECTLY TO THE SKIN)	10 x 15 cm (4" x 6")
715	COLD THERAPY COMPRESS	15 x 26 cm (6" x 10")

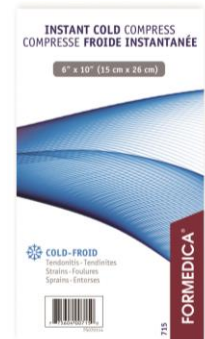

Instructions for use

1. Fold top to bottom to "POP" inner fluid bag
2. Shake
3. Apply to injured area

WARNING :

- Use only as directed. Do not apply directly to skin. Always wrap in towel or cloth before application.
- Keep away from eyes when activating the cold pack. Should contents become exposed to skin or eyes, flood with water immediately and consult a physician.
- Do not apply to open wound. Always test temperature for comfort level before application.
- For external use only. The unattended use of a cold pack by children, elderly or incapacitated person is dangerous.
- Harmful if swallowed. If accidentally swallowed, drink several glasses of water, call your local Poison Control Center immediately.
- Discard if punctured or leaking.
- Keep out of reach of children.